
Operating Instructions | Instrucciones de funcionamiento

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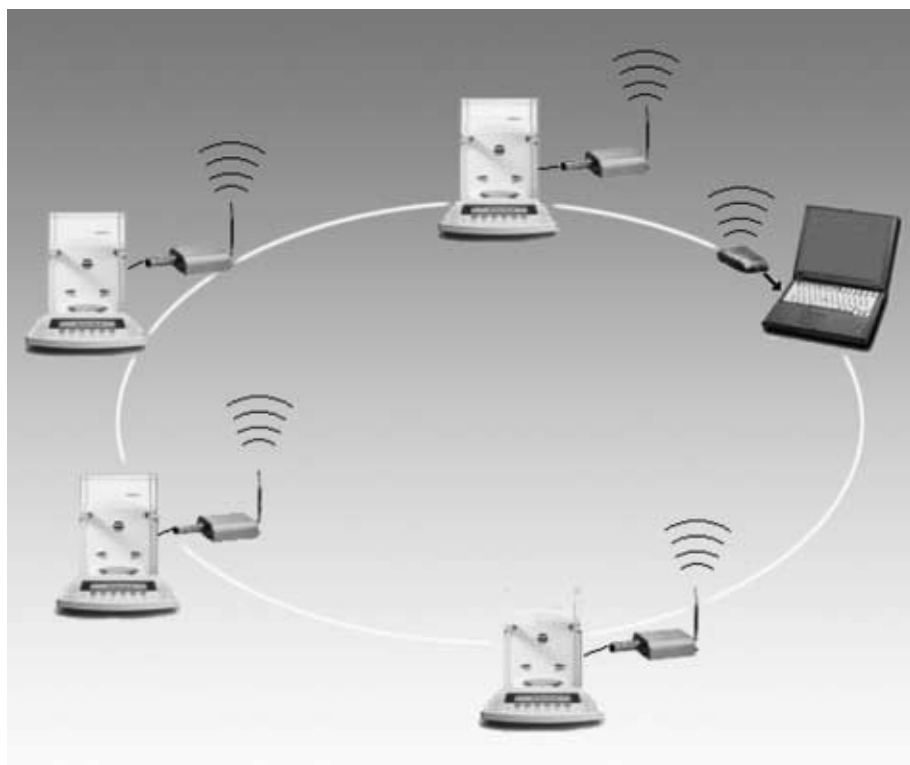English – page 3

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Contents

- 3 **Intended Use**
- 4 **Warnings and Safety Precautions**
- 7 **Getting Started**
- 11 **Installing the YBT01**
- 15 **Installing the YD001C-BT (Option B2)**
- 17 **Establishing Data Connections between the YBT01 or YD001C-BT (Option B2) and the F8T001v2 or F8T012 from Belkin Corporation**
- 35 **Care and Maintenance**
- 37 **Specifications**
- 38 **Certificates**

The following symbols are used in these instructions:

- indicates required steps
- indicates steps required only under certain conditions
- > describes what happens after you have performed a certain step
- indicates an item in a list
- ⚠ indicates a hazard

Make sure you observe the warnings and safety precautions in their entirety during installation and operation, as well as while performing maintenance and repair work on the equipment. It is important that all personnel using the equipment understand this information, and have access to it at all times. The warnings and safety precautions can be supplemented by the equipment operator. Make sure all operating personnel are informed of any additions to these instructions.

Intended Use

The wireless interface from Sartorius enables wireless data communication, using **Bluetooth®** wireless technology, between a measuring instrument and a remote device. This remote device may be a PC, for example, equipped with a **Bluetooth®** USB adapter. With the wireless interface, you no longer need data cables to connect your RS-232 data interfaces.

The YBT01 is a **Bluetooth®** RS-232 adapter that can be connected to the RS-232 port on a weighing instrument.

YD001C-BT transmitter (Option B2) is also a **Bluetooth®** RS-232 adapter, for installation on the UniCOM interface in a Combics indicator.

The YBT01 and YD001C-BT (Option B2) interfaces have a transmit power of up to 100 milliwatts (in accordance with **Bluetooth®** wireless technology specification V1.1/Class 1).

Once the radio connection has been established, measured values can be processed on the PC using, for example:

- WinScale®, or
- SPEXXIS® from Sartorius.

Warnings and Safety Precautions

Legal Notices

- **Bluetooth®** is a registered trademark of the **Bluetooth®** Special Interest Group (**Bluetooth®** SIG, Inc.). The **Bluetooth®** word mark and logos are owned by the **Bluetooth®** SIG, Inc. and any use of such marks by Sartorius is under license. Other trademarks and trade names are those of their respective owners.
- The **Bluetooth®** USB adapter model F8T001v2 and **Bluetooth®** USB adapter F8T012 are products of the Belkin Corporation. All rights reserved. The F8T001v2 and F8T012 **Bluetooth®** USB adapters are available from Sartorius AG under the model name “YBT02.” Sartorius reserves the right to integrate adapters from other manufacturers, however, in place of the F8T001v2 or F8T012 **Bluetooth®** USB adapter.
- Windows® 98, Windows® 2000 and Windows® XP are registered trademarks of Microsoft Inc.
- WinScale® and SPEXXIS® are registered software products of Sartorius AG.

Copyright

The programs and products named in these instructions are subject to proprietary, copyright and other laws. The software products from the various manufacturers may not be changed, reverse engineered, or modified through assimilation. The purchaser/end user does not acquire any of the rights named above at any time.

The software may not be reproduced or copied. The purchaser/end user may operate the software for his/her own purposes only, and may not permit use of the software by third parties, neither free of charge nor in return for payment.

If you have any problems with the software supplied, please inform the software vendor. Sartorius AG reserves the right to supply updated software with this product. No liability is accepted for software installed prior to the purchase of this product.

The manufacturer of the **Bluetooth®** USB adapter is responsible for the content of the compact disk (CD) supplied with the **Bluetooth®** USB adapter, which includes the operating instructions, the programs and software, legal notices and installation instructions.

Safety Information

- Please read these operating instructions carefully before using your wireless interface.
The manufacturer is not liable for damage resulting from failure to observe these warnings and safety precautions
 - Read the operating instructions supplied with the Bluetooth® USB adapters (on the enclosed CD) carefully as well, and make sure you observe all warnings and safety precautions therein.
 - △ The wireless interface must not be used in hazardous areas/locations in which explosive gases or dusts are present, nor in areas potentially endangered by explosive material.
 - △ The wireless interfaces may not be operated in legal metrology, nor in locations in the medical industry
 - △ When the wireless interface is in use, the risk of disturbance caused by other approved devices that use the same frequency cannot be ruled out.
If you need assistance, contact your Sartorius dealer or the Sartorius Service Center.
 - Have the equipment inspected by a trained technician at appropriate intervals for correct functioning and safety. Make sure the applicable accident prevention regulations are observed by all operating personnel.
 - △ Always make sure the equipment is disconnected from AC power before performing any installation, cleaning, maintenance or repair work. If the equipment housing is opened by anyone other than persons authorized by Sartorius, this will result in forfeiture of all claims under the manufacturer's warranty. Use only original Sartorius spare parts.
 - △ Disconnect the equipment from power before connecting or disconnecting peripheral devices.
 - △ The casing on all connecting cables, as well as the casing on wires inside the equipment housing, is made of PVC. Chemicals that corrode these materials must be kept away from these cables.
 - △ Operating the equipment beyond the specifications stated in the Declaration of Conformity and in these instructions is not permitted, and constitutes use of the equipment for other than its intended purpose. All specifications stated in the Declaration of Conformity and in these operating instructions must be strictly observed. The device meets the EC Directives listed in the Declaration of Conformity. Any installation work that does not conform to the instructions in this manual will result in forfeiture of all claims under the manufacturer's warranty.
 - △ Do not expose the equipment to aggressive chemical vapors or to extreme temperatures, moisture, shocks, or vibration. The allowable operating temperature range during operation is 0°C to +40°C (+32°F to +104°F).
- CE !**
- △ The YBT01 and YD001C-BT (Option B2) interfaces are based on Bluetooth® wireless technology. This wireless technology operates in a frequency range from 2.4 to 2.4835GHz, which is not standardized through Europe. Always be sure to observe the regulations applicable at the place of installation when using this equipment.
 - △ In France, the wireless interface may be used only within buildings or enclosed spaces.

Installation

△ Proceed with extreme caution when using pre-wired RS-232 connecting cables, as the pin assignments may not be compatible with Sartorius equipment.

- The operator shall be solely responsible for any damage or injuries that occur when using cables not supplied by Sartorius.

△ If there is visible damage to any of the equipment or power cords, disconnect the equipment from power and lock it in a secure place to ensure that it cannot be used for the time being.

IP Rating:

△ The IP65 protection rating for the YBT01 is ensured only if installation is performed correctly in accordance with the instructions and acknowledged technological standards. The IP67 protection rating for the YD001C-BT (Option B2) is ensured only when it is installed in the Combics CIS indicator correctly in accordance with the instructions and acknowledged technological standards.

The YD001C-BT must be installed and tested by a certified technician.

Using the Equipment in the EU:

- The model YBT01 and YD001C-BT (Option B2) equipment is intended for use only in the following countries:

Austria, Belgium, Denmark, Finland, France (indoor use only), Germany, Greece, Iceland, Ireland, Liechtenstein, Luxembourg, Norway, Portugal, Spain, Sweden, Switzerland, The Netherlands, United Kingdom.


Getting Started

Storage and Shipping Conditions

Do not expose the YBT01 or YD001C-BT (Option B2) to extreme temperatures, shocks, vibration or moisture.


Unpacking the Equipment

- After unpacking the equipment, check it immediately for any visible sign of damage resulting from rough handling during shipment.
- If any sign of damage is visible, refer to the section entitled “Safety Inspection” in the chapter entitled “Care and Maintenance.” It is a good idea to save the box and all parts of the packaging until you have successfully installed your wireless interface. Only the original packaging provides the best protection for shipment. Before shipping the equipment, disconnect all cables and unscrew the antenna from the YBT01 or YD001C-BT (Option 2) interface to prevent damage during shipment.

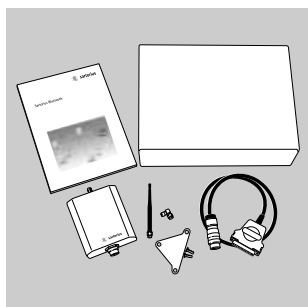
Equipment Supplied for YBT01

- YBT01 (RS-232 adapter)
- Connecting cable
- Angle plug
- Antenna
- Bracket for wall mounting
- Screws (2) for installing wall-mounting bracket
- Packaging
- Plastic bags (2)
- Operating instructions (this manual)

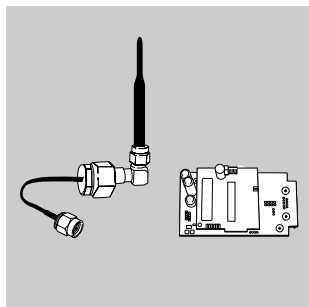Equipment Supplied for YDO01C-BT (Option B2)

- YDO01C-BT (RS-232 adapter)
- Connecting cable
- Angle plug
- Antenna
- Seal
- Nut (1)
- 2 printed circuit boards (PCBs)
- Packaging (not with factory-installed Option B2)
- Plastic bag (not with factory-installed Option B2)
- Equipment designation labels (2) (not with factory-installed Option B2)
- Operating instructions (this manual)

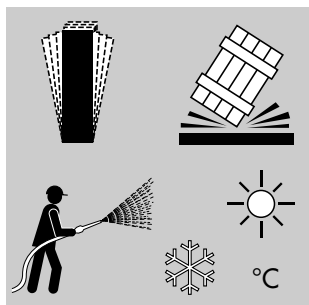Installation Instructions

The YBT01 and YDO01C-BT (Option B2) wireless interfaces are designed to provide reliable results under normal ambient conditions. When choosing locations to set up your equipment, please observe the following to ensure fast and precise operation of the YBT01 or YDO01C-BT (Option B2) wireless interface:

- Position the interfaces for optimum signal transmission between devices.

Important Note:

Solid obstacles such as walls can restrict the transmission range.

- Avoid placing the equipment in close proximity to a heater or otherwise exposing it to heat or direct sunlight.
- Avoid locations that are subject to excessive vibration
- Protect the wireless interfaces from highly concentrated aggressive chemical vapors.
- Avoid locations subject to excessive moisture
- The wireless interfaces must not be used in hazardous areas/locations, nor in areas in which explosive gases or dusts are present.
- The equipment may not be operated in areas potentially endangered by explosive material, nor in locations in the medical industry.

Software from Sartorius

Sartorius offers two software programs for fast and easy processing of the data and control commands transmitted between connected devices:

- WinScale
 - SPEXXIS
- Each of these programs comes with complete instructions for simple installation on the PC used with your wireless interfaces.

Configuring the Weighing Instrument

Data transmission parameters in the RS-232 data interface on the YBT01 are configured as follows:

- Transmission speed (baud rate): 9600 bit/s
- Parity: Odd
- Stop bits: 1
- Data bits: 7

Data transmission parameters in the RS-232 data interface on the YD001C-BT (Option B2) are configured as follows:

- Transmission speed (baud rate): 9600 bit/s
- Parity: Odd
- Stop bits: 1
- Data bits: 8

Thus the data interface on the weighing instrument must be adapted for data transmission to the YBT01 or YD001C-BT (Option B2) wireless interface. Please observe the following when configuring the weighing instrument:

The weighing instrument models suitable for wireless data communication use either the SBI (Sartorius Balance Interface) or the XBPI (Xtended Balance Protocol Interface) data transfer protocol. Which protocol is used in your weighing instrument is indicated in the

operating instructions for the particular instrument.

The XBPI protocol generally operates with the following parameters:

- Transmission speed (baud rate): 9600 bit/s
- Parity: Odd
- Stop bits: 1
- Data bits: 8

With the SBI protocol, the transmission speed, parity and number of stop bits can be configured. The following Sartorius instruments also permit configuration of the number of data bits:

- GENIUS (ME models)
- MA100 and MA50 (electronic moisture analyzers)
- Combics (CISL... models) (indicators and complete scales)
- Combics (CIS... models) (indicators and complete scales)

Refer to the installation and operating instructions for the particular instrument for details on which parameters can be configured.

To enable communication with the YBT01 or YD001C-BT (Option B2) wireless interface, the data interface on the weighing instrument must be configured accordingly.

To use the SBI protocol and a YBT01 with its factory-set parameters, configure the weighing instrument as follows:

- Transmission speed (baud rate): 9600 bit/s
- Parity: Odd
- Stop bits: 1
- Data bits: 7

Note:

The model YBT01 **Bluetooth®** RS-232 adapter with factory-set interface parameters is not suitable for data communication using the XBPI protocol. To use the weighing instrument and the YBT01 with XBPI, or with SBI in a configuration other than the factory-set configuration of the serial interface, the desired parameters must be configured in the serial interface of the YBT01. The serial interface of the **Bluetooth®** RS-232 model YBT01 adapter can be configured at the factory in accordance with specified parameters, or by one of our service technicians in accordance with customer specifications.

To use the SBI protocol and a YD001C-BT (Option B2) with factory-set parameters, configure the weighing instrument as follows:

SBI parameters:

- Transmission speed
(baud rate): 9600 bit/s
- Parity: Odd
- Stop bits: 1
- Data bits: 8

To use the XBPI protocol and a YD001C-BT (Option B2) with factory-set parameters, configure the weighing instrument as follows:

- Select "XBPI" in the scale or balance operating menu

Note:

To use the weighing instrument and the YD001C-BT with SBI in a configuration other than the factory-set configuration of the serial interface, the desired parameters must be configured in the serial interface of the YD001C-BT (Option B2). The serial interface of the

Bluetooth® model YD001C-BT adapter can be configured at the factory in accordance with specified parameters, or by one of our service technicians in accordance with customer specifications.

The model YBT01 **Bluetooth®** RS-232 adapter is designed for data transfer with Sartorius premium laboratory balances and analyzers of the following series:

- GENIUS (ME models)
- Masterpro (LA models),
- MA100 und MA50 (electronic moisture analyzers)
- Expert (LE models)

The model YBT01 **Bluetooth®** RS-232 adapter is designed for data transfer with Sartorius industry scales of the following series:

- Combics (CISL... models) (indicators and complete scales)

The model YD001C-BT **Bluetooth®** RS-232 adapter (Option B2) is designed for data transfer with Sartorius industry scales of the following series:


- Combics (CIS... models) (indicators and complete scales)

Installing the YBT01


If you wish to install the YBT01 on the wall, attach the enclosed bracket for wall mounting before proceeding with the following instructions (see “Installing the Bracket for Wall Mounting,” below).

Installing the YBT01


- Connect the threaded angle plug (SMA adapter) to the jack on the back of the YBT01 transmitter

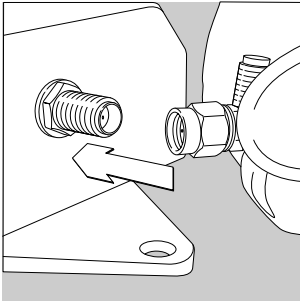
- Use an open-ended wrench (8 mm) to tighten the plug carefully

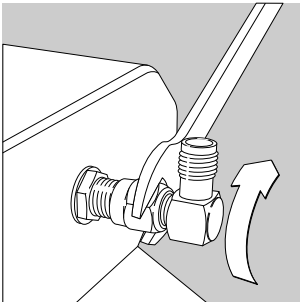
- Plug the antenna into the jack on the angle plug and tighten it

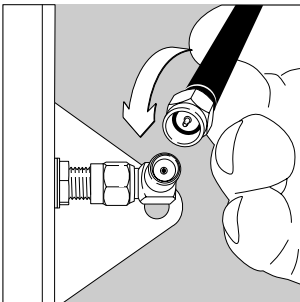
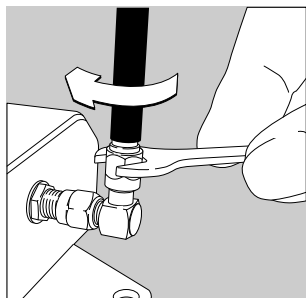
- Use an open-ended wrench (8 mm) to tighten the antenna carefully

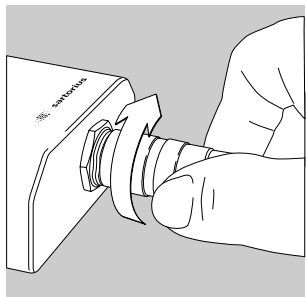
- Connect the cable (RS-232) to the threaded jack on the front of the YBT01

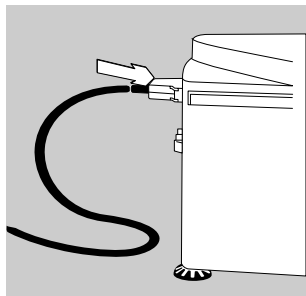- Connect the other end of the cable to the RS-232 data interface on the weighing instrument and tighten the screws on the connector

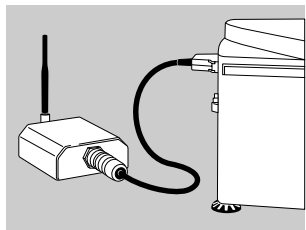
- This completes installation of the transmitter; the YBT01 is ready to operate

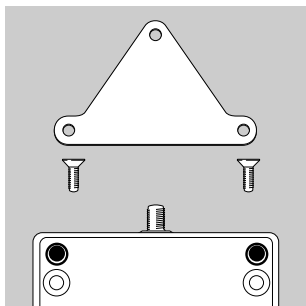
Installing the Bracket for Wall Mounting


- Remove the bracket and screws from the packaging

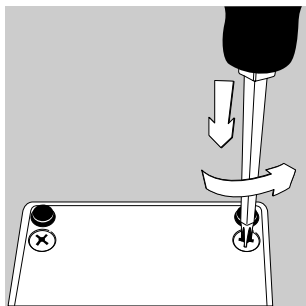
- Remove the two screws nearest the antenna jack on the YBT01. Save these two screws; e.g., in the plastic bag supplied.

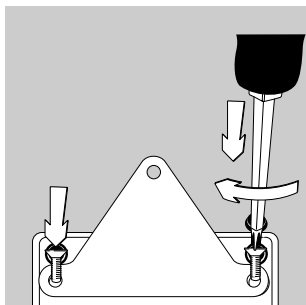
- Use the screws supplied with the bracket to attach the bracket to the YBT01 (see illustration)

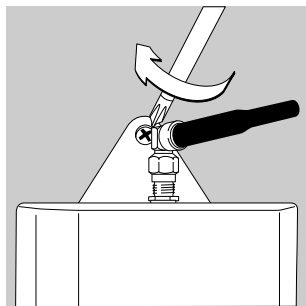- Mount the YBT01 on the wall

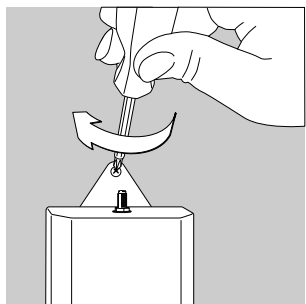
- Alternatively, you can install the antenna without the angle plug. To do this, connect the antenna directly to the jack on the YBT01.

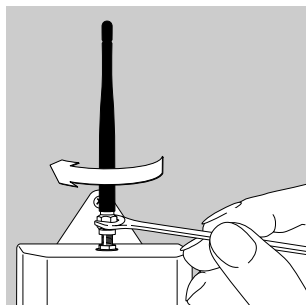
- Use an open-ended wrench (8 mm) to tighten the antenna carefully. Store the angle plug in a safe place.


Installing the YD001C-BT (Option B2)

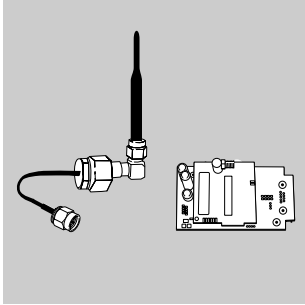
- ⚠ Because the installation of the YD001C-BT (Option B2) entails opening the indicator housing, it must be performed by a technician with relevant training (i.e., a service technician). Failure to observe the applicable restrictions and regulations when opening the indicator housing will result in forfeiture of all claims under the manufacturer's warranty.

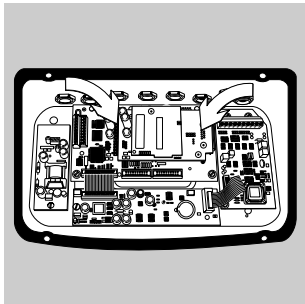
- Attach the PCBs

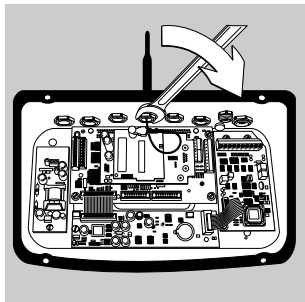
- Remove the protective cap from the back of the indicator housing
- Pull the antenna cable through the opening
- Fasten the nut by hand and then use an open-ended wrench (19 mm) to tighten the antenna connector carefully
- Plug the antenna into the angle plug and tighten it by hand
- Use an open-ended wrench (8 mm) to tighten the antenna carefully

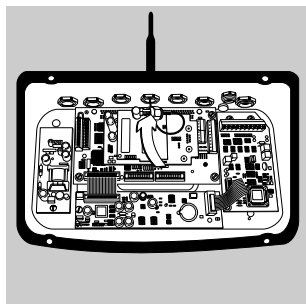
- Attach the antenna cable to the threaded connector on the YDO01C-BT adapter

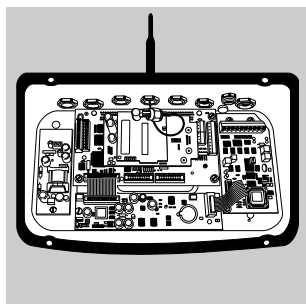- Close the indicator housing

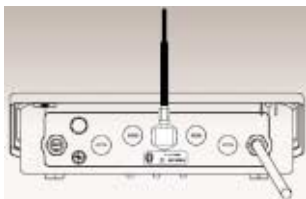- Affix the enclosed ID label to the rear panel of the indicator as shown:

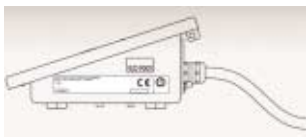Affix the warning label to the side panel of the indicator next to the CE marking:

Establishing Data Connections between the YBT01 or YD001C-BT (Option B2) and the F8T001v2 or F8T012 from Belkin Corporation

Note:

Up to four connections can be used simultaneously.

My Bluetooth Places

Installing F8T001v2 or F8T012

- Install F8T001v2 or F8T012 in accordance with the installation instructions provided on the CD.

Note:

A detailed installation manual for the F8T001v2 is included on the enclosed CD.

- Once installation has been completed, double-click on “My Bluetooth Places.”
 - > The “Initial Bluetooth Configuration Wizard” opens.

Note:

The Bluetooth Configuration Wizard starts automatically only the first time you open “My Bluetooth Places.”


- Add or delete checkmarks in the boxes to indicate which shortcut(s) to Bluetooth functions you wish to install.
- Click on [Next] to continue.
 - > The “General” dialog opens.

- Enter a name for the computer; e.g., the name already assigned by the network
- Select the computer type from the list; e.g., “Laptop”
- Click on [Next] to continue.
 - > The “Initial Bluetooth Configuration Wizard” screen is shown again, with a different message than shown on the first screen.


- Click on [Next] to continue.
 - > The “Bluetooth Service Selection” page opens.

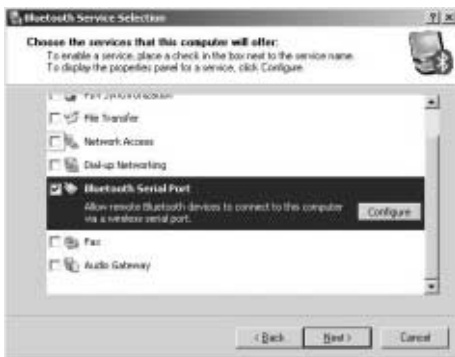- Add or delete checkmarks in the boxes to select the required Bluetooth services.

Note:

Only the “Bluetooth Serial Port” service is required for a data connection between the YBT01 or YD001C-BT and the F8T001v2 or F8T012. Activate the corresponding checkbox to select this service.
- Click on the [Configure] button.
 - > Windows installs the software automatically.
 - Afterwards, the “Bluetooth Properties” page opens.

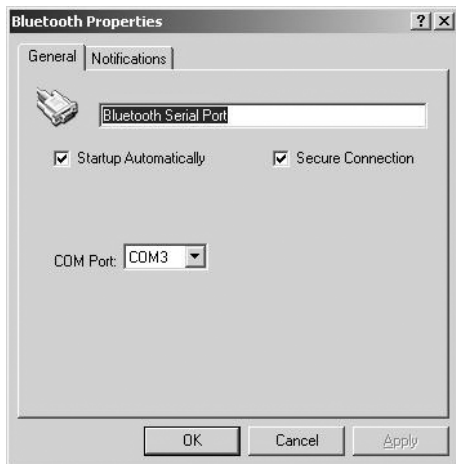- Edit the name of the **Bluetooth** port as needed.
- Activate the “Startup Automatically” option to have the connection made on this port automatically.

Note:

Depending on the receiving device, automatic start might not be possible.

- Select the “Secure Connection” option to require input of a PIN for authentication of the two devices before the connection can be made.

Note:

The factory-set (default) PIN for the YBT01 and the YD001C-BT is 0000. You can have the Sartorius Service department change this number, if desired.

- Click on the “Notifications” tab.
- > This opens the Notifications page.

- Activate/deactivate the checkboxes and edit the input field as desired to define the manner of notification given when a connection is established.

Note:

Click on the [Browse] button to search for a different audio file. Click on the [▶] button to play the audio file.

- Click on [OK] to confirm your settings.
 - > The “Bluetooth Service Selection” screen opens again.

- Click on [Next] to continue.
 - > The screen shown here is displayed

- Click on the [Skip] button

- Click on the [Finish] button to conclude the initial configuration of your Bluetooth environment.
- > The basic Bluetooth® configuration is now completed.

Establishing a Connection

Configuring the Virtual Serial COM Ports

Note:

Configuration of multiple virtual COM ports is required only if you wish to use more than one connection at a time. Once you have configured the COM ports your settings are saved; you do not need to configure the ports again the next time you open “My Bluetooth Places” or re-start the computer.

My Bluetooth Places


- Double-click on “My Bluetooth Places.”
 - > The “My Bluetooth Places” window opens.


- Open the “Bluetooth” menu.


- Select “Advanced Configuration” from the “Bluetooth” menu.
 - > The “Bluetooth Configuration” panel opens.

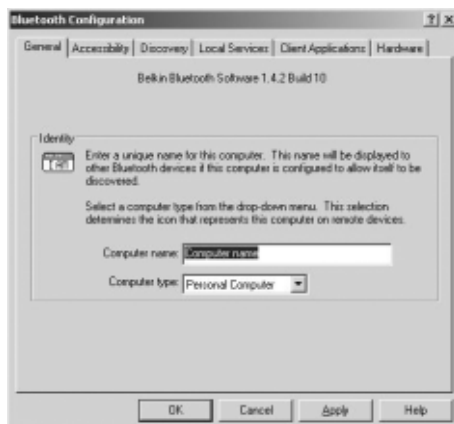
- Click on the “Client Applications” tab.

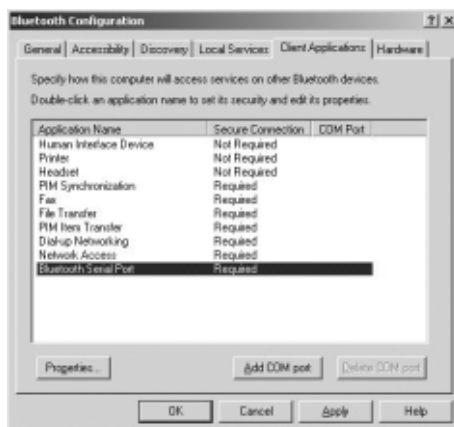- Click on “Bluetooth Serial Port” to select it.
 - Click on the [Properties] button.
- > Windows installs a virtual COM port for the application automatically. The name of the COM port is assigned automatically (in this example, COM4). The “Bluetooth Properties” page opens.

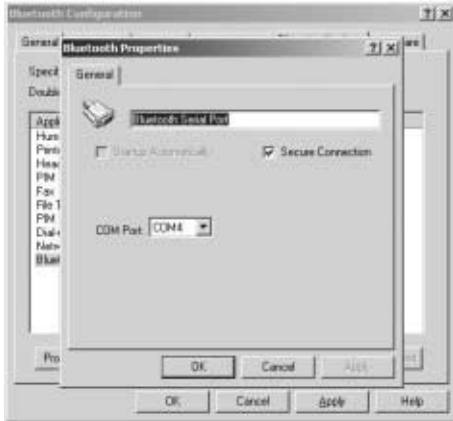
- Edit the name of the Bluetooth port in the input field, if desired.
 - Select the “Secure Connection” option to require input of a PIN for authentication before a connection can be made on this port.
 - Click on the [OK] button to confirm the settings for the selected COM port.
- > The “Bluetooth Properties” page closes.

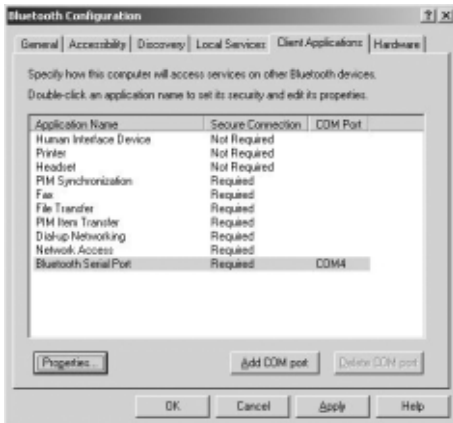Creating Additional COM Ports

- Click on the [Add COM port] button.
- > Windows automatically installs an additional virtual COM port for the application. The name of the COM port is assigned automatically (in this example, COM5). The “Bluetooth Properties” page opens.

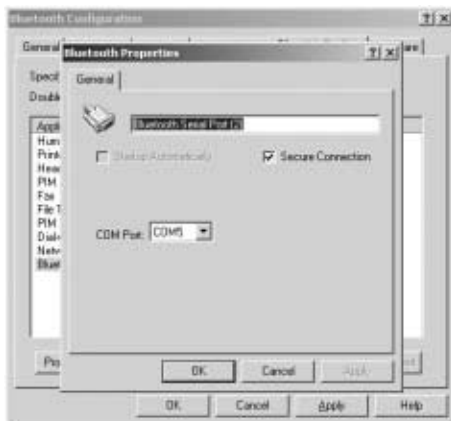- Edit the name of the **Bluetooth** port in the input field, if desired.
- Select the “Secure Connection” option to require input of a PIN for authentication before a connection can be made on this port.
- Click on the [OK] button to confirm the settings for the selected COM port.
- > The “**Bluetooth Properties**” page closes.

- Repeat these steps to configure up to four separate serial COM ports for **Bluetooth** connections.
- Note:**
These ports can be used by applications such as “WinScale” or “SPEXXIS” to access up to four weighing instruments.
- Click on the [Apply] button to save your **Bluetooth** configurations.
- > The computer is configured with your settings.
- Click on [OK] to close the “**Bluetooth Configuration**” panel.

Connecting to other Bluetooth Devices

Search for and identify devices within range of the Bluetooth USB adapter.

- Open the “Bluetooth” menu.
- Select the “Search for Devices” item from the “Bluetooth” menu.
- > The program searches for devices that are equipped with Bluetooth® wireless technology and are within range of the F8T001v2 or F8T012
- > The devices detected are displayed with a symbol and a device name in the window pane on the right.

Note:

If more than one device is detected with the same name, the devices cannot be unambiguously identified from the device names. For clear identification, you need to check the Device Properties for the device address and compare this number with the ID numbers printed on each device. The ID number of the YBT01 is printed on the manufacturer's ID label, with is affixed to the bottom of the device. The ID number of the YD001C-BT (Option B) is printed on the manufacturer's ID label affixed below the antenna connector (see "Installing the YD001C-BT (Option B2)" in this manual for details).

- Select the device you wish to identify (by clicking on it) and then press the mouse button on the right ("right-click"); this opens a device-specific menu ("popup menu").
- Select the "Properties" menu item.
 - > The "<Device Name> Properties" window opens.
 - > The identification number is shown here as the "Device Address" in hexadecimal code.
- Click on [OK] to close the "<Device> Properties" window.

Searching for Available Services on the Remote Device


- Select the remote device and then right-click on it to open the popup menu.
- Select the “Discover Available Services” menu item.
- > The services available on the remote device are listed. This is also known as the device profile.

Note:

Depending on the computer's operating system, the symbols for devices and services may differ from those shown here. Devices that are equipped with **Bluetooth®** wireless technology can be connected to one another only if they support the same profiles.

- Select the “serial port” service and then right-click on it to open the popup menu for services.

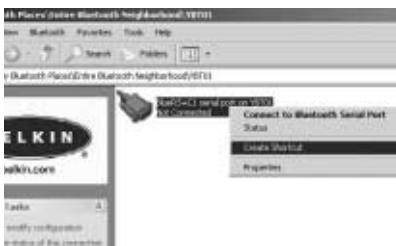- Select the “Create Shortcut” menu item.
- > A message is displayed indicating that a shortcut has been created at the highest folder level of the “My Bluetooth Places” window.

Note:


The next time you start “My Bluetooth Places” you can use this shortcut to establish the connection, without having to search for the remote device and detect its services. To establish a connection using the shortcut, make sure the remote device is switched on and is within range of the Bluetooth USB adapter.

- Click on [OK] to close the message window.
- Click on the [Back] button.
- You can use the same procedure to create shortcuts for all of the remote devices with which you wish to establish data connections.
- Click on “My Bluetooth Places” in the “Other Places” section on the left,

or:

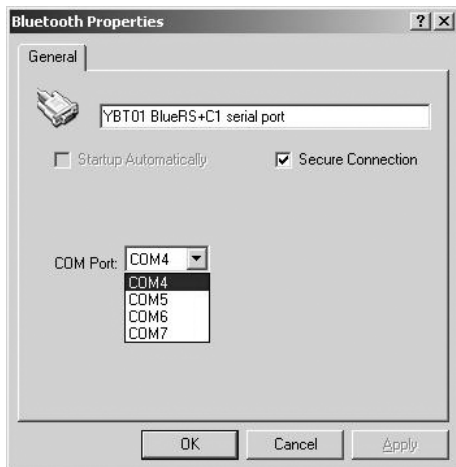- Edit the name of the shortcut in the input field, if desired. Select the “Secure Connection” option to require input of a PIN for authentication before a connection can be made using this shortcut.

- Select the desired COM port from the list.

Note:

The remote device can communicate with application software over the COM port you select here. The COM ports must be set up prior to opening the “Properties” page.


- Click on [Apply] and then [OK] to confirm the settings for the selected shortcut.

> The “Bluetooth Properties” page closes.

- Repeat this procedure to assign COM ports to the other connection shortcuts.

Note:

If a particular COM port is assigned to more than one shortcut, the connections started by these shortcuts cannot be used simultaneously. For connections that you wish to use simultaneously, separate COM ports must be assigned to the respective shortcuts. You can use up to four YBT01 and YDO01C-BT (Option B2) connections simultaneously.

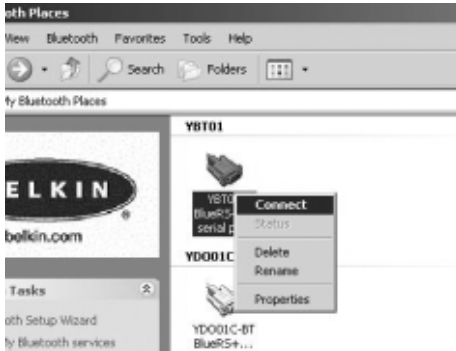Activating a Data Connection

- Select the connection shortcut and then right-click on it to open the popup menu for the shortcut.

- Select the “Connect” menu item.

- > A status message is displayed, as well as a message indicating that PIN code must be entered.


Note:

The PIN input prompt is displayed only if the “Secure Connection” option is selected in the shortcut properties.

- Click on the **Bluetooth®** logo in the system tray of the task bar.

- > A dialog box for PIN input opens.

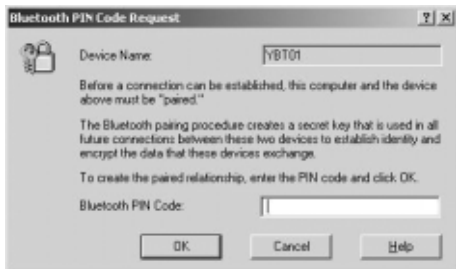- Enter the PIN in the input field and click on [OK] to confirm.

Note:

The factory-set (default) PIN for the YBT01 and the YD001C-BT (Option B2) is 0000. You can have the Sartorius Service department change this number, if desired.

- > A message window opens. The message provides information on the configuration of the active port. The application software (e.g., WinScale or SPEXXIS) can use the established connection on the port indicated here.

Click on the [Properties] button to open the “Bluetooth Properties” window. The properties cannot be edited while the connection is active.

Click on the [Disconnect] button to close down the connection.

- Click on the [Close] button to close the “Bluetooth Connection Status” window.

Note:

For more information, refer to the functions offered in the “Help” menu in the “My Bluetooth Places” window.

Care and Maintenance

Service

Regular servicing by a Sartorius technician will ensure continuous correct functioning of your wireless interface. Sartorius offers service contracts with regular maintenance intervals ranging from 1 month to 2 years. The maintenance interval depends on operating conditions and tolerance requirements.

Repairs

- ⚠ Disconnect defective equipment from power immediately. Repair work must be performed by authorized Sartorius service technicians using original spare parts. Repairs performed by untrained persons may result in considerable hazards for the user. If the equipment is still under warranty, send the entire indicator to the factory for repairs.

Safety Inspection

Safe operation of the equipment is no longer ensured when:

- there is visible damage to the device or cable,
 - the device no longer functions properly,
 - the device has been stored for a relatively long period under unfavorable conditions (e.g., excessive moisture), or
 - the equipment has been subjected to rough handling during shipment.
- If there is any indication that safe operation of the device is no longer warranted:
- disconnect the equipment from power immediately.
 - Lock the device in a secure place to ensure that it cannot be used for the time being, and
 - notify your nearest Sartorius Service Center or the International Technical Support Unit based in Goettingen, Germany.

Cleaning

- ⚠ Handle the YBT01 and YD001C-BT/Option B2 wireless interfaces with care.
- ⚠ Make sure no liquid enters the equipment housing.
- ⚠ Do not use aggressive cleaning agents (solvents or similar agents).
- The accessible parts (antenna, adapter, housing) can be cleaned by wiping with a cloth slightly moistened with a mild detergent.
- ⚠ Using jets of water or compressed air to clean the wireless interfaces is not permitted.

Cleaning Stainless Steel Surfaces

Clean all stainless steel parts regularly. Wipe down stainless steel parts gently (e.g., the YBT01 housing) using a damp cloth or sponge. You can use any household cleaning agent that is suitable for use on stainless steel. Then wipe the equipment thoroughly to remove all residues. Afterwards, allow the equipment to dry. If desired, you can apply oil to the cleaned surfaces as additional protection.

Note:

Do not apply detergents or cleaning agents to ID labels or printed adhesive tags!

Corrosive Environment

- Remove all traces of corrosive substances on a regular basis.

Recycling

The packaging protects the equipment from damage during transport. After successful installation of the equipment, you should return this packaging for recycling. For information on recycling options, including recycling of used equipment, contact your municipal waste disposal center or local recycling depot.

Serial Number Coding

The month and year of manufacture are encoded in the serial number as follows:

Y M M x x x x x

Y	Year
1	2000–2006
2	2007–2013
3	2014–2020
4	2021–2027
5	2028–2034
6	2035–2041
7	2042–2048
8	2049–2055
9	2056–2062

The first digit represents a 7-year period (see the table above); the next 2 digits represent the month. Months are numbered consecutively over the entire 7-year period, starting at 13. Thus the number representing the month also indicates the specific year of manufacture.

2000 13–24

2001 25–26

...

Example:

113xxxxx -> January 2000

The individual devices are numbered consecutively in the last 5 digits, starting from 00000 again at the beginning of each month.

Specifications

YBT01

Transmit power	In acc. with Class 1
Specification	Bluetooth® V1.1
Profiles supported	Serial port profile
Data transmission speed	1200 to 115,200 bit/s
Operating temperature range	0°C to +40°C (+32°F to +104°F)
Type of IP protection	IP65
Dimensions (D × W × H)	121 mm × 84 mm × 32 mm (without connecting cable, antenna and bracket for wall mounting)

YD001C-BT (Option B2)

Transmit power	In acc. with Class 1
Specification	Bluetooth® V1.1
Profiles supported	Serial port profile
Data transmission speed	1200 to 115,200 bit/s
Operating temperature range	0°C to +40°C (+32°F to +104°F)

CE Declaration of Conformity to Council Directives 1999/5/EC, 89/336/EEC and Recommendation 1999/519/EC

**The telecommunications terminal equipment
YBT01 and YDO01C-BT**

meets the requirements of the test standards listed below when used in accordance with regulations and in conjunction with the associated power supplies, auxiliary peripheral devices and installation equipment listed in Annex A2 (for a list of the individual type designations and technical descriptions, see Annex A1).

1. Essential requirements of Article 3 of Council Directive 1999/5/EC

1.1 Sources of 1999/5/EC: Official Journal of the European Communities No. 2002/C62/03

- EN 60950 Safety of information technology equipment including electrical business equipment
- EN 301 489-1 Electromagnetic compatibility and radio spectrum matters (ERM); electromagnetic compatibility (EMC); standard for radio equipment and services
Part 1: Common technical requirements
- EN 301 489-17 Part 17: Specific conditions for 2,4 GHz wideband transmission systems and 5 GHz high performance WLAN equipment
- EN 300 328-2 Electromagnetic compatibility and radio spectrum matters (ERM); wideband transmission systems; data transmission equipment operating in the 2,4 GHz ISM band and using spread spectrum modulation techniques
Part 2: Harmonized European Standard (EN) covering essential requirements under article 3.2 of the R&TTE Directive

2. Electromagnetic compatibility


2.1 Source of 89/336/EEC: Official Journal of the European Communities No. 2002/C62/02

- EN 61326-1 Electrical equipment for measurement, control and laboratory use / EMC requirements
Part 1: General requirements
Limitation of emissions: Residential areas, Class B
Defined immunity to interference: Minimum test requirements, continuous unmonitored operation

3. Protection against exposure to electromagnetic fields

- EN 50371 Generic standard to demonstrate the compliance of low power electronic and electrical apparatus with the basic restrictions related to human exposure to electromagnetic fields (10 MHz - 300 GHz) - General public

Sartorius AG
37070 Goettingen, Germany
2003

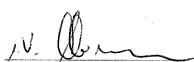W. Obermann
(Senior Vice President, R&D
Electronics Engineering
Mechatronics Division)

Dr. K. Klein
(Senior Vice President
Mechanical Engineering
Mechatronics Division)

Bluetooth® Qualification Body (BQB)

Michael Klos
CETECOM ICT Services GmbH, Untertürkheimer Str. 6-10, 66117 Saarbrücken, Germany

QPN Qualified Product Notice

Number of annexes: 1

Reference number of the QPN: **QPN-T813816P-01-BQB**

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Product name: **Bluetooth Interface E/I**

Product ID(s): **YBT01, YDO01C-BT, Option B2**

Product category: **Computer Accessory**

Product type: **PROD**

Hardware version: **1.0**

Software version: **BTA1000, V1.124**

Short Product Description: **Bluetooth Adapter external / internal, Class 1**

Profiles supported: **GAP, Serial-DevA, Serial-DevB, SDAP-LocDev**

This Qualified Product Notice (QPN) is only valid in conjunction with Annex 1 and the Compliant Portion Declaration (CPD).

The Qualified Product Notice is an essential part of the *Bluetooth* qualification program and shows the acceptance of this product as a qualified *Bluetooth* product. On the basis of the reviewed submitted compliance folder for the above mentioned product we certify that the product has passed all requirements to be listed as a qualified *Bluetooth* product in accordance to chapter "6.3 Re-Marketing of Listed Products" of the *Bluetooth* Qualification Program Reference Document (QPRD) v1.0 and the specifications for *Bluetooth* wireless technology as described in annex 1 and in the Compliant Portion Declaration (see separate document on the QPL with the reference number CPD-T813816P-01-BQB).

Saarbrücken,
Place, Date of Issue 2003-11-10
D-66117 Saarbrücken

Signed by Michael Klos, *Bluetooth* Qualification Body - BQB –
CETECOM ICT Services GmbH
Phone: +49 (0) 681 598-3722
Fax: +49 (0) 681 598-3775
E-mail: Michael.Klos@ict.cetecom.de

The *Bluetooth* word mark and logos are owned by the Bluetooth SIG, Inc. and any use of such marks by CETECOM is under license. Other trademarks and trade names are those of their respective owners.

Bluetooth® Qualification Body (BQB)

Michael Klos

CETECOM ICT Services GmbH, Untertürkheimer Str. 6-10, 66117 Saarbrücken, Germany

CPD

Compliant Portion Declaration


Number of annexes : 1

Reference number of the CPD: **CPD-T813816P-01-BQB**
Reference number of the QPN: **QPN-T813816P-01-BQB**
Assessment Date: **2003-01-31**
Listing Date: **2003-11-10**
Applicant details: **Sartorius AG**
Weender Landstraße 94-108
37075 Göttingen, Germany
Manufacturer/ member: **See above**
Product contact person: **Swen Weitemeier**
e-mail: swen.weitemeier@sartorius.com
Product name: **Bluetooth Interface E/I**
Product ID(s): **YBT01, YDO01C-BT, Option B2**
Product category: **Computer Accessory**
Product type: **PROD**
Hardware version: **1.0**
Software version: **BTA1000, V1.124**
Short Product Description: **Bluetooth Adapter external / internal, Class 1**
Profiles supported: **GAP, Serial-DevA, Serial-DevB, SDAP-LocDev**

This Compliant Portion Declaration (CPD) is only valid in conjunction with Annex 1 and the Qualified Product Notice (QPN) with the reference number QPN-T813816P-01-BQB.

On the basis of the reviewed submitted compliance folder for the above mentioned product we certify that the mandatory and optional features of the Bluetooth Protocols and Profiles as described in Annex 1 and as implemented in the above mentioned Product are in compliance to the Bluetooth Qualification requirements.

Saarbrücken,
Place, Date of Issue: 2003-11-10
D-66117 Saarbrücken

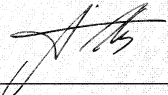Signed by Michael Klos, Bluetooth Qualification Body - BQB -
CETECOM ICT Services GmbH
Phone: +49 (0) 681 598-8722
Fax: +49 (0) 681 598-8775
E-mail: Michael.Klos@ict.cetecom.de

The Bluetooth word mark and logos are owned by the Bluetooth SIG, Inc. and any use of such marks by CETECOM is under license. Other trademarks and trade names are those of their respective owners.

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Saarbrücken,
Place, Date of Issue 2003-11-10
D-66117 Saarbrücken

Signed by Michael Klos, *Bluetooth* Qualification Body - BQB –
CETECOM ICT Services GmbH
Phone: +49 (0) 681 598-8722
Fax: +49 (0) 681 598-8775
E-mail: Michael.Klos@ict.cetecom.de

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
On the basis of the reviewed submitted compliance folder for the above mentioned product we certify that the mandatory and optional features of the *Bluetooth* Protocols and Profiles as described in Annex 1 and as implemented in the above mentioned Product are in compliance to the *Bluetooth* Qualification requirements.

Saarbrücken,
Place, Date of Issue: 2003-11-10
D-66117 Saarbrücken

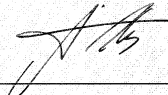Signed by Michael Klos, *Bluetooth* Qualification Body - BQB –
CETECOM ICT Services GmbH
Phone: +49 (0) 681 598-8722
Fax: +49 (0) 681 598-8775
E-mail: Michael.Klos@ict.cetecom.de

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Copyright by Sartorius AG,
Goettingen, Germany.
All rights reserved. No part
of this publication may
be reprinted or translated in
any form or by any means
without the prior written
permission of Sartorius AG.
The status of the information,
specifications and illustrations
in this manual is indicated
by the date given below.
Sartorius AG reserves the
right to make changes to the
technology, features,
specifications and design of the
equipment without notice.

Status:
February 2006, Sartorius AG,
Goettingen, Germany